

CORIO CD-1201F Refrigerated / heating circulator

Refrigerated / heating circulators of the CORIO CD series distinguish themselves with a great price-to-performance ratio. Durable and dependable, they are suitable for a wide range of standard and routine work in both laboratories and industrial environments.

The CORIO CD-1201F refrigerated circulator is exceptionally versatile with a bath opening of 35×41 cm (width x depth) and a bath depth of 30 cm. The resulting bath volume of 50 liters makes it possible to control the temperature of large objects, or many smaller objects, directly in the bath.


This cooling machine uses natural, environmentally-friendly refrigerant and was developed with a focus on energy efficiency. This means significant savings of up to 70% with respect to the operating costs for numerous applications, which also means faster amortization of the procurement cost. At the same time, the lower energy consumption positively contributes to climate protection.



Product features

- Models for internal and external applications
- Bright, white, easy to read display
- Very quiet
- USB connection
- Integrated pump connection M16×1
- · easy switching between internal and external circulation
- · Bath lid and drain tap included
- removable ventilation grid
- Refrigeration unit without side vents
- Class III (FL) according to DIN 12876-1
- Space-saving cooling coil design yields more usable space in the bath tank

Pump capacity

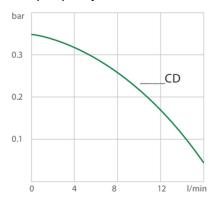

Medium: Water



Performance values

230V/50Hz (Schuko Plug - CEE 7/4 Plug Type F)			
Heating capacity kW	2		
Viscosity max. cSt	50		
Pump capacity flow pressure I/min	15		
Pump capacity pressure bar	0.35		
Power consumption A	11		

Order number including voltage version			9012718.N1.0	9012718.N1.03			
Cooling capacit	y (Ethanol)						
°C	20	10	0	-10	-20	-30	-40
kW ¹	1.25	1.2	1.1	0.85	0.54	0.3	0.05
Refrigerant stag	ge 1						
Refrigerant		R1270					
Filling weight g		90					
Global Warming R1270	g Potential for	2					
Carbon dioxide	equivalent t	0.00018					

¹ Performance specifications measured in accordance with DIN 12876. Cooling capacities up to 20 °C measured with ethanol; over 20 °C with thermal oil unless otherwise specified. Performance specifications apply at an ambient temperature of 20 °C. Performance values may differ with other bath fluids.

Technical data

Available voltage versions		Cooling		
Order No.	9012718	Cooling of compressor	1-stage Air	
Available voltage versions:				
9012718.N1.03	230V/50Hz (Schuko Plug - CEE 7/4 Plug Type F) (R1270)			
9012718.N1.03.chn	230V/50Hz (CN Plug) (R1270)			
9012718.N1.04	230V/50Hz (UK Plug Type BS1363A) (R1270)			
9012718.N1.05	230V/50Hz (CH Plug Type SEV 1011) (R1270)			
9012718.N1.13	230V/60Hz (Nema N6-20 Plug) (R1270)			
9012718.N1.01	100V/50-60Hz (Nema N5-20 Plug) (R1270)			
9012718.N1.02	115V/60Hz (Nema N5-20 Plug) (R1270)			

Bath		Other	
Bath tank	Stainless steel	Classification	Classification III (FL)
Bath cover	integrated	IP Code	IP 21
Usable bath opening cm (W x L / D) 35 x 40 / 30		Pump function	Pressure Pump
		Pump type	Immersion Pump
Electronics		Dimensions and volumes	
Interfaces	USB	Weight kg	68



Temperature control	PID1	Barbed fittings inner diameter mm	8/12 mm
Absolute temperature calibration	1 Point Calibration	Total dimensions cm (W \times L \times H)	45 x 64 x 95
Temperature display	LED	Filling volume I	48 56
Temperature setting	Keypad	Pump connections	M16x1 male

Temperature values	
Setting the resolution of the temperature display °C	0.1
Working temperature range °C	-40 + 100
Temperature stability °C	±0.03
Ambient temperature °C	+5 +40
Temperature display resolution °C	0.01 0.1

Included in delivery

2 Barbed fittings for tubing 8 and 12 mm ID. (Pump connections M16x1 male) $\,$

All Benefits



Brilliant.

Very bright display makes it easy to read even from a distance.



Switch on. And off you go.

Intelligent operating concept. Ready for operation with just a few quick and easy steps.



Timer. Integrated.

CORIO circulators include an integrated timer function. When the set time has elapsed, a signal sounds and the device switches off. Setting range: 0 ... 999 minutes.



Everything at the front.

All operating controls and safety functions are accessed easily and comfortably from the front.



Exact.

You can rely on it. PID1 control and 'Active Cooling Control' make the new CORIO precise and perfect.



Locked in.

The lockable power plug guarantees safe connection. More process safety.



Connection. Easy.

Inclined pump connections (M16×1) facilitate the connection of applications. Each unit includes 2 barbed fittings of 8/12 mm diameter each.



Early warning system for low liquid level

Maximum safety for applications, optical and audible alarm, allows user to refill bath fluid before the unit shuts down



Connectivity.

Remote control made easy. CORIO CD circulators feature a USB connection.